

Science. Results. Trust.

"elōs" technology has been phenomenal in my practice. They combine energies in such a way that we can provide exceptionally effective treatments in the safest possible way with the greatest level of comfort and without downtime, and that's what our patients value about elōs®."

Terry Jacobson, M.D., Washington, USA

Technical Specifications

Applicator	Spectrum	Light Fluence	RF Energy	Treated Area	Pulse Rep Rate	Cooling
Sublative	N/A	N/A	Up to 62 mJ/pin	Standard 64: 12 x 12 mm Focal 44: 11 x 3 mm	1 Hz	N/A
Sublime O	700-2000 nm	Up to 6 W/cm ²	Up to 120 J/cm ³	12 x 8 mm	1 Hz	10°C
SR O	580-980 nm	Up to 45 J/cm ²	Up to 25 J/cm ³	25 x 12 mm	1 Hz	5°C
SRA O	470-980 nm	Up to 45 J/cm ²	Up to 25 J/cm ³	25 x 12 mm	1 Hz	5°C
AC O	400-980 nm	Up to 18 J/cm ²	Up to 25 J/cm ³	25 x 12 mm	1 Hz	5°C
LV	915 nm	Up to 140 J/cm ²	Up to 100 J/cm ³	8 x 5 mm	1 Hz	5°C
LVA	915 nm	100-350 J/cm ²	Up to 100 J/cm ³	8 x 2 mm	1 Hz	5°C
Motif LHR-LHR Mode	810 nm	Up to 50 J/cm ²	Up to 50 J/cm ³	15 x 12 mm	3 Hz	5°C
Motif HR - HR Mode O	680-980 nm	Up to 45 J/cm ²	Up to 25 J/cm ³	25 x 12 mm	1 Hz	5°C
Motif IR - IR Mode O	915 nm	Up to 70 J/cm ²	Up to 100 J/cm ³	8 x 5 mm	1 Hz	5°C

Motif Modes	Frequency	Total Energy Per Pulse	Total Energy Per Second
Motif LHR - Motif Mode	10 Hz	14 J	140 J
Motif HR - Motif Mode	3 Hz	16 J	48 J
Motif IR - Motif Mode	10 Hz	2.5 J	25 J

Platform Weight	87 kg / 132 lbs	ative
Platform Platform Size W x D x H]	54 x 52.5 x 140 cm / 21.5 x 22 x 55.2 in	dica Cre
Electrical Specifications	100-240 VAC	⊚ Me







Dubai Medica Trading LLC Health Care City Office:

Tel: + 971 4 43 70 669

Abu Dhabi Medica Trading LLC


Tel: + 971 2 65 81 619 Tel: + 966 1 46 54 113

Saudi Arabia Medica / RIC

Lebanon Medica ME s.a.l.

Tel: +961 471 0973 e-mail: info@medica.bz e-mail: info@medica.bz e-mail: info@medica.bz e-mail: info@medica.bz Built for the Next Generation Multi Platform



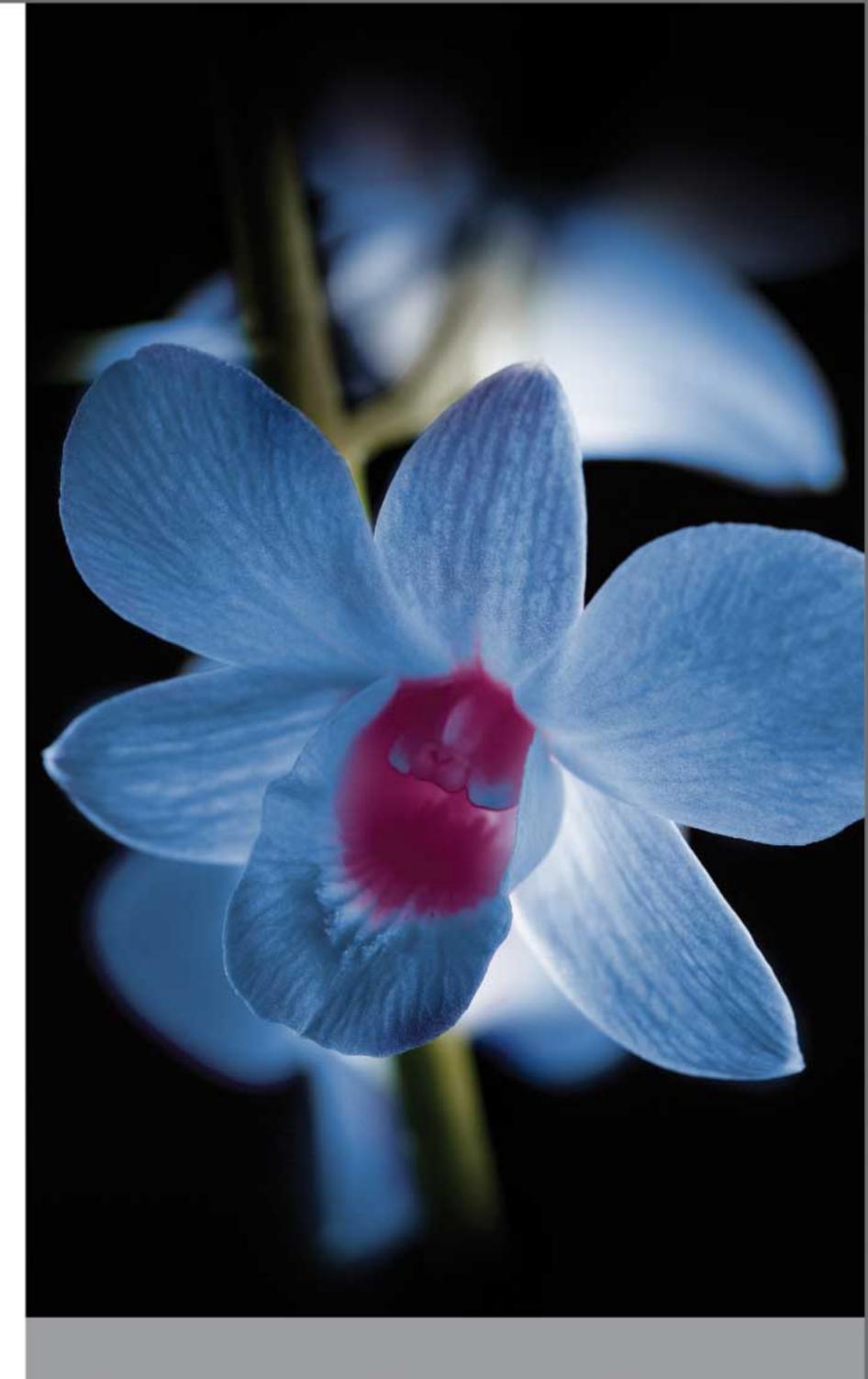

Science. Results. Trust.



$OOSPUS^{m}$

Offers more than 8 in-demand applicators for more than 30 treatments that will maximize your practice revenue.

LASER Applicators

[MOTIF LHR]

Pain Free Hair Removal 810 nm



Advanced Vascular Lesions Treatment 915 nm



MOTIF IR

Focused Wrinkle Reduction 915 nm



SUBLATIVE

Bipolar Fractionated Radio Frequency



SUBLATIVE RF

Fractional Skin Treatment Increase in Collagen and elastin with minimal epidermal distraction

IPL (Intense Pulsed Light) Applicators



[SUBLIME]

Facial Contouring & Laxity Treatment 700-2000 nm



[MOTIF HR]

Pain Free Hair Removal 680-980 nm



SR/SRA_ (skin rejuvenation)

Advanced Vascular & Pigmented Lesions 470-980 nm




eMotif™

eStyle™Motif™





Complete Facial Rejuvenation

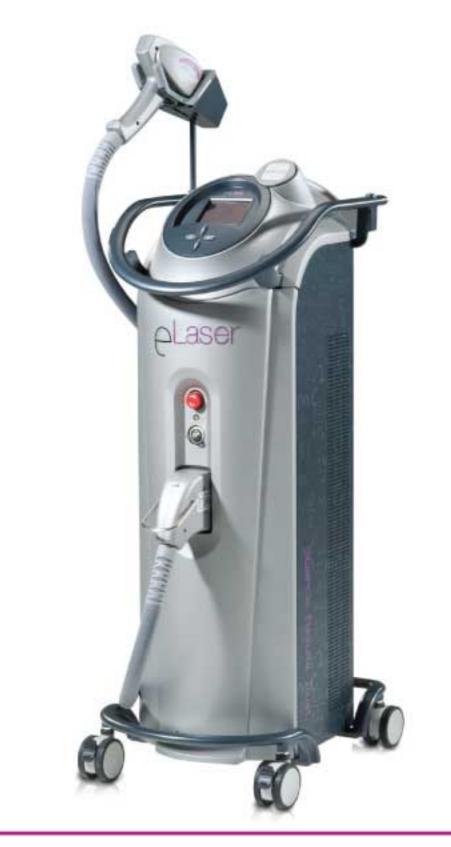


Complete Laser + RF Platform



Complete IPL + RF Platform



Fast.
Safe.
Pain Free.
Effective.
Combination

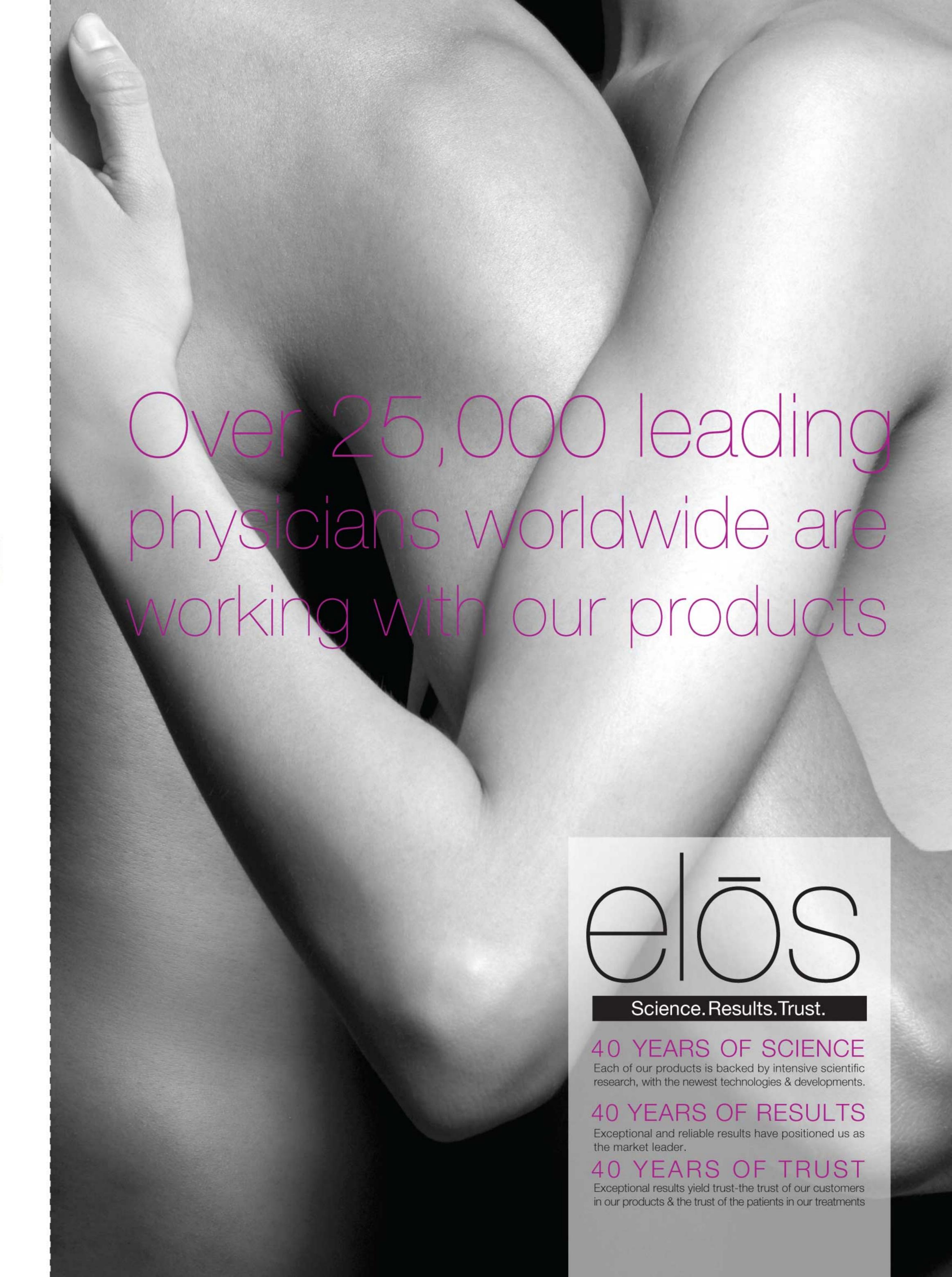



Finally, a multi-application system for professionals!

The new elos Plus™ Complete Aesthetic Work Station utilizes

- elōs Technology a combination of optical energy and bi-polar radio frequency (RF) for safe and effective treatments on all skin types.
- Upgradeable, with 8 different treatment applicators
- Intuitive user interface with easy to use guided treatment options
- Motif Mode™ for fast, effective and "Pain Free" hair removal and skin treatments
- Sublative™ fractionated bi-polar radio frequency (RF) for increased collagen with minimal epidermal disruption
- triniti™ a total skin rejuvenation protocol in the same visit



Features and Benefits

elōs Plus™

	Intuitive Touch Screen Interface	Improved ease of use
System	Guided Treatment Modes	Quick Reference Guide is embedded in the Interactive screen
	Plug & Play software upgrades	Easy to upgrade applications
	Cooling capacity increased versus e-max	System works faster

	Combination of Sublime & Sublative	Multi-level treatment of the skin for complete facial rejuvenation in one visit
	Integration of IFS (Intelligent Feedback System)	Ensures precise RF energy delivery regardless of skin impedance variations
Clinical	The Triniti treatment protocol	SR/SRA: Color correction Sublime: Facial contouring & laxity treatment Motif IR: Focused wrinkle treatment
	Full Sublative iD support on the Sublative RF applicator	Maximizes treatment options and results
	Addition of Motif™ modes to 3 applicators	Enables "Pain Free" hair removal and skin treatments
	High output frequency of up to 3 Hz in LHR (normal) mode & 10 Hz in Motif mode	Fast treatment times
	New RF electronics measure RF more accurately	Improves bad coupling and partial contact issues

Safety & Patient Comfort		
	Active Dermal Monitoring™	ISM / ISL Monitors the skin impedance changes throughout each pulse & cuts pulse when reaching set limit Strong safety profile
	Integrated sapphire contact cooling	Increased patient comfort and added safety
	Short / long pulse - Enables selection sub-pulse frequency	Added safety and patient comfort
	Optimized tip geometry on Motif HR applicator	Increased patient comfort

Industrial Design	Single turn applicator connectors	Easy to connect / disconnect applicators
	Water purity sensor	Informs user if water needs to be changed
	Water bottle located in the front of the system	Easy access to water bottle
	Tilting screen	Various viewing angles



Provides total facial renewal with no downtime!

SR/SRA

Skin Rejuvenation / Acne, Hyperpigmentation, Superficial Veins

- Treatment of pigmented and vascular lesions.
- Overall improvement in skin appearance.

SR Applicator





Refore

Photos courtesy of Ronald C. Russo, MD

Sublime

Skin Tightening and Contouring

- Non-invasive laxity treatment.
- Uses the elos combination of gentle pulsed infrared light and bi-polar radio frequency (RF)
- Stimulates long term collagen synthesis, resulting in an immediate & long term firmed appearance.

Sublime Applicator





Before

Post 2 Treatments

Photos courtesy of Ruthie Amir, MD

Motif IR

Focused Wrinkle Treatment

- The strongest and most focused treatment of deep wrinkles, skin texture and pore size.
- Safer & more effective treatments on all skin types.
- "Pain Free" treatment of superficial skin texture and pore size.

Motif IR





Before

Post 1 Treatment

Photos courtesy of John Shieh, MD



- ☑ Complete Full Range Applications
- ✓ Safe & Effective on All Skin & Hair Types
- ✓ Innovative Feedback System for Optimal Results!

Applicators / Before & After Results!


Motif LHR / HR

Fast and Pain Free Hair Removal

- Achieves permanent hair removal of pigmented and minimally pigmented hair on all skin types, including tanned skin.
- Enables a high energy hair removal mode, as well as "Pain Free" Motif mode.
- Provides highly selective targeting of the hair follicle and shaft while protecting the surrounding tissue.

Motif Hair Removal



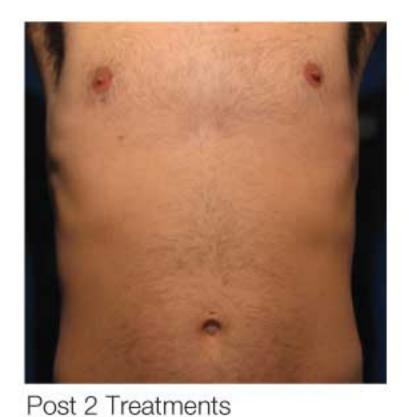

Before

Photos courtesy of Ruthie Amir, MD

Sublative Rejuvenation

Increase in Collagen

- Produces a significant increase in both collagen and elastin with minimal epidermal disruption.
- Initiates a strong wound healing response in the dermal matrix while also minimizing the damage to the epidermis.
- Assures a smoother, rejuvenated appearance and a reduction in fine lines, wrinkles, thin skin laxity, skin irregularities and striae.

Sublative





Before

Post 2 Treatments

Photos courtesy of Stephen Bassett, MD

LV / LVA

Leg Veins (Superficial and Deep)

- Easily treats vascular lesions in 1-3 treatments.
- Treats vascular lesions and leg veins up to 4mm.
- Perfect for advanced vascular and smaller resistant vessels including fine telangiectasias.
- Provides controlled selective thermolysis to vascular lesions, telangiectasias and leg veins.
- Delivers maximum targeted energy to the treatment area.

LV Applicator





Before


Post 3 Treatments

Photos courtesy of Cyrus Chess, MD


The following applicators are available with the system:

- Sublative Rejuvenation / Increase in Collagen
- Motif LHR & Motif HR / Fast & Pain Free Hair Removal
- SR / SRA / Skin Rejuvenation, Acne, Hyperpigmentation, Superficial Veins
- Sublime / Skin Tightening and Contouring
- Motif IR / Focused Wrinkle Treatment
- LV / LVA / Leg Veins (Superficial and Deep)



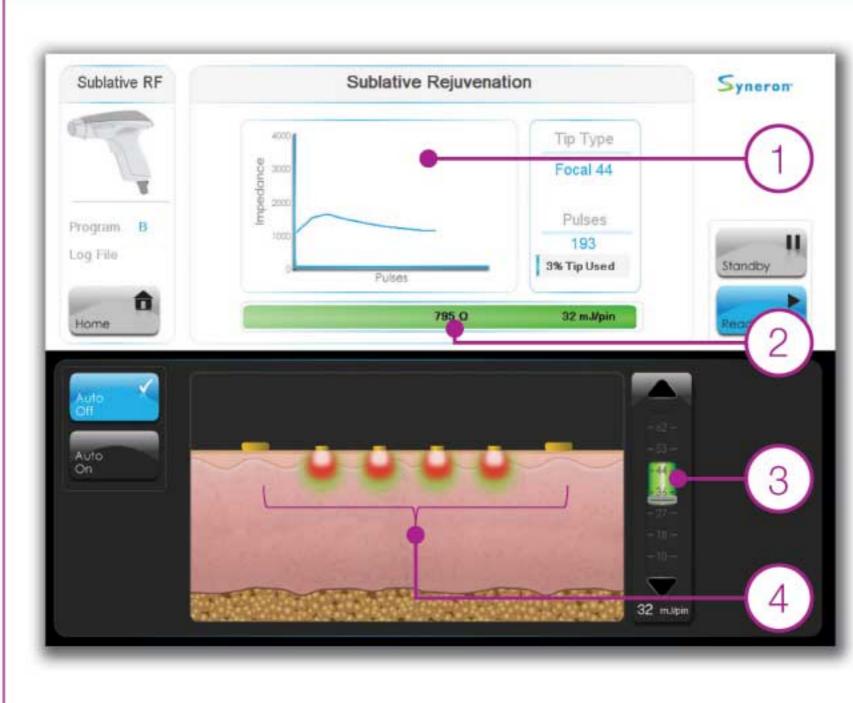

- 1- Intelligent Feedback Graph (IFG) Shows a moving average of pulse to pulse impedance.
- 2- Skin Feedback Meter (SFM) Shows skin condition and energy delivered.
- 3- Full Energy selection bar With highlighted program range.
- 4- Scaled injury pattern Shows depth and width of sublative thermal damage.

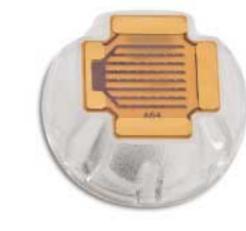

Sublative iD Standard 64 pins 100 & 400 pulses per tip



Sublative iD Focal 44 pins 200 pulses per tip

The Sublime & Sublative Advantage

Re-texture and reverse the signs of aging skin with the innovative combined treatment of Sublime and Sublative applicators. Comprehensive facial rejuvenation is achieved through multi-level enhancement of collagen.

Combined Treatment





Photos courtesy of César Arroyo Romo, MD


Combined Treatment






Photos courtesy of Ruthie Amir, MD

Science of Sublime

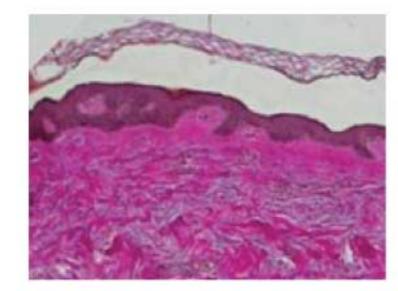

Pre Treatment

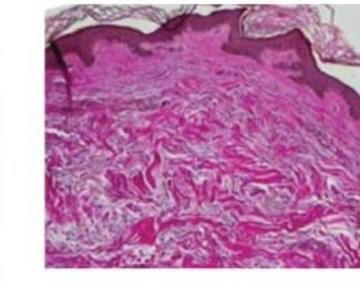



Deep dermal remodeling 3 months post treatment


Su-blime [suh-blahym], noun.

- 1. Non-invasive wrinkle treatment
- 2. elōs technology for safe and effective treatments
- 3. Fast treatment, immediate visible results

Science of Sublative

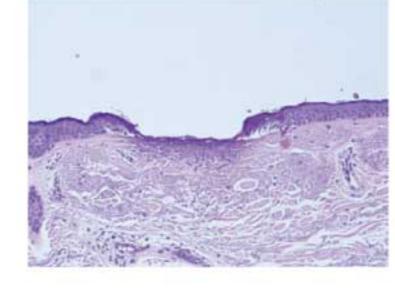

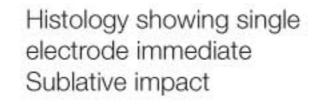

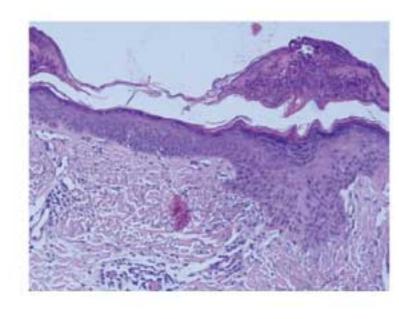

Histology showing rapid healing and strong dermal remodeling only 4 days post a single

Sub-la-tive [sah-bla-tiv], noun.

- 1. A new category of fractional skin treatments
- 2. High dermal impact, low epidermal disruption for all skin types
- 3. Minimal downtime, proven efficacy, satisfied patients



- ✓ Complete Full Range Applications
- Safe & Effective on All Skin & Hair Types
- ✓ Innovative Feedback System for Optimal Results!