

4 Technologies  
in 1 PLATFORM

**ADVANCED  
BODY  
SCULPTING  
PROCEDURES\***

- RF | Radio Frequency
- DD | DeepTDepression
- IR | InfraRed
- MR | Mechanical Rollers



	<i>VelaShape II</i>	
	Body Applicator	Contour Applicator
RF Power	up to 60W	up to 23W
Infra Red Power	up to 35W	up to 20W
Treatment Area	40x40 mm	30x30 mm
Light Spectrum	700 - 2000 nm	
Vacuum	Pulsed	
Platform Weight	27 Kg	
Platform Size	36 x 55 x 80 cm	

\*For professional use only



ISO 9001:2008  
Certificate No. 34-8-07-QMS 0039

**m** **medica**

[www.medica.bz](http://www.medica.bz)

**Dubai  
Medica Trading LLC**  
Head Office:

Tel: + 971 4 28 28 338  
e-mail: [info@medica.bz](mailto:info@medica.bz)

**Dubai  
Medica Trading LLC**  
Health Care City Office:

Tel: + 971 4 43 70 669  
e-mail: [info@medica.bz](mailto:info@medica.bz)

**Singapore  
Medica Pte Ltd**

Tel: + 65 6733 9689  
e-mail: [info-sg@medica.bz](mailto:info-sg@medica.bz)

**Saudi Arabia  
Medica / R/C**

Tel: + 966 1 46 54 113  
e-mail: [info@medica.bz](mailto:info@medica.bz)

**Lebanon  
Medica ME s.a.l.**

Tel: + 961 4 71 09 73  
e-mail: [info-me@medica.bz](mailto:info-me@medica.bz)

**NEW!**

vela-series

# Vela™

Powered by elōs™ Technology



*VelaShape II™*

**MEDICAL  
AESTHETIC SOLUTIONS**

A revolutionary approach for medical-aesthetic treatment solutions

The first, class II\* medical device to receive:

**FDA clearance for  
CIRCUMFERENTIAL REDUCTION  
BODY CONTOURING  
& CELLULITE REDUCTION**

**CE marked for  
BODY CONTOURING  
& CELLULITE REDUCTION**

# Synergy of 4 technologies in 1 PLATFORM

elōs™ technology

**Class II\* system:**

A Class II laser/light system is safe because the blink reflex of the human eye will limit the exposure to no more than 0.25 seconds. Which means that there will be no possibility of eye damage under any circumstances.



**VISIBLE RESULTS  
in as Few as 4 TREATMENTS**

- 4 Treatments protocol for visible results
- Treat abdomen, buttocks or thighs in less than 20 minutes
- Proprietary elōs™ technology for optimal results
- More power for shorter treatment sessions

**The EXTRAORDINARY Science  
of Velashape II™ Synergistic Use of Energies**

The Velashape II™ device employs four treatment modalities:

- Infrared light (IR) heats the tissue up to 5 mm depth
- Radio Frequency (RF) heats tissue from 2 to 20 mm depth
- Vacuum provides improved circulation and a precise demarcated heating effect
- Mechanical massage provides lymphatic drainage (dermal and hypodermal)

**Biological Impact of Treatment**

Increased metabolism leading to the desired anabolic and catabolic alterations at the predetermined skin components and depths.


**Mechanical MANIPULATION**

- Stimulates lymphatic drainage and blood perfusion
- Facilitates fibroblast activity
- Reduces viscosity of fat cell clusters
- Promotes vasodilatation and extravasations of oxygen and nutrients
- Enables heating at different depths

**HEATING**

- Enhances blood perfusion and increases oxygen dissociation from oxy-hemoglobin
- Facilitates fibroblast activity
  - Increases metabolism of fat cells
  - Improves skin texture

## How elōs™ Works...

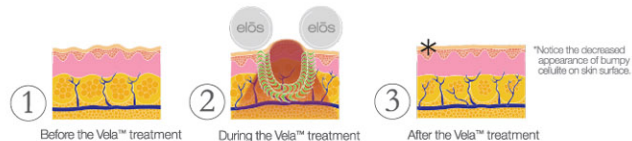

### "Body Applicator"

■ LARGE APPLICATOR & DISPOSABLE HEAD



### "Contour Applicator"

■ SMALL APPLICATOR & DISPOSABLE HEAD



- The Vela™ System is enhanced with 2 applicators:
  - The **Body applicator** is used for large areas such as the buttocks, thighs and abdomen.
  - The **Contour applicator** is used for small areas such as the arms, flanks and double chin.

## Before & After

